



# AioCare Spirometry Better Diagnostics Better Lives

A completely mobile and flexible approach



## Hardware

Lung diagnostics  
& environment  
monitoring



## Software

Mobile app &  
secure cloud-  
based data  
platform



Time efficient spirometry and  
diagnosis for experts and general  
practitioners alike.



Wireless and truly mobile.



Holistic patient treatment  
with smart home monitoring.

# AioCare Spirometry Components

## Full Spirometry


AioCare registers inspiratory & expiratory flow volume curves corrected by built-in real-time BTPS & displays key parameters, as well as SpO<sub>2</sub>/HR.

## Calibration free

Be ready for spirometry at any time.

## Simple Solution

Seamless IT integration. Wireless connection & charging.

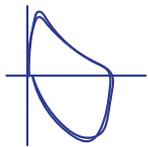

### Data

Quality & reproducibility assured through AI algorithms.



### Hygiene

Protect patients, staff & device from cross-contamination




MicroGard™ II filters provide 99.999% viral & bacterial efficiency against cross contamination (Nelson Test Report 1003754).<sup>1</sup>



# AioCare Fully connected Smart home monitoring

## Doctor app

Result instantly available for review and guidance

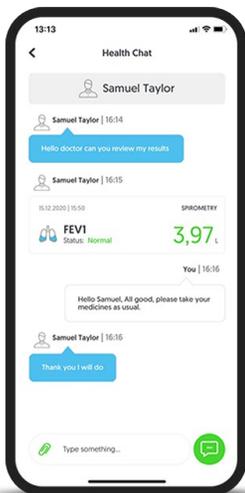


AIOCARE SERVER /  
DOCTOR PANEL  
AMAZON AWS  
ISO 27001



## Patient app

Ability to chat with health care professionals

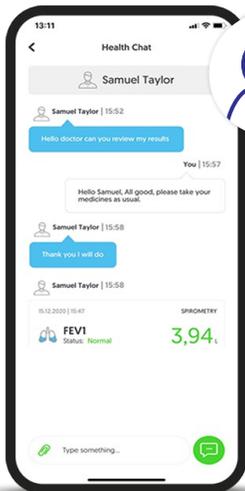

Contact initiated by patient. Measurement results sent



Answer sent by doctor



# Why AioCare?



## Confirmed quality

A cutting edge, class IIa medical device (MDD), certified by the notified body TUV Nord. Validated by an observational study in a sample of 80 patients.



## Accessibility

A small, portable device that can be taken anywhere. Safety tested by UL. Easy to use and inexpensive. Intuitive interface for creating patient records automatically saves test results, enables statistical analysis and management of the disease history.



## 4 easy steps to your first results

High accuracy of test results identical to those conducted with the medical equipment available at medical centres.



## Warranty and support

2-year warranty and customer support line.

## AioCare. Better lives.

AioCare means higher quality care for your patients. You give your patients a device to self-monitor and take care of their own health in a safe way, under your care.

Patients know their data is automatically shared with their doctor. It saves hours of telephone consultations, urgent visits and avoidable emergency interventions.

## Stay connected with your patients regardless of their geographical location.


An electronic journal enables the patient to avoid filling out tables by hand. The results of your AioCare patients are instantly visible in your panel. Patients are guided through the correct test technique through clear instructions on the patient App to ensure spirometry measurements are performed in a technically acceptable fashion. AioCare supports the user with additional sensors and a virtual guide to perform spirometry measurements in a technically correct way.

# AioCare Spirometry

## A completely mobile and flexible approach



### AioCare App

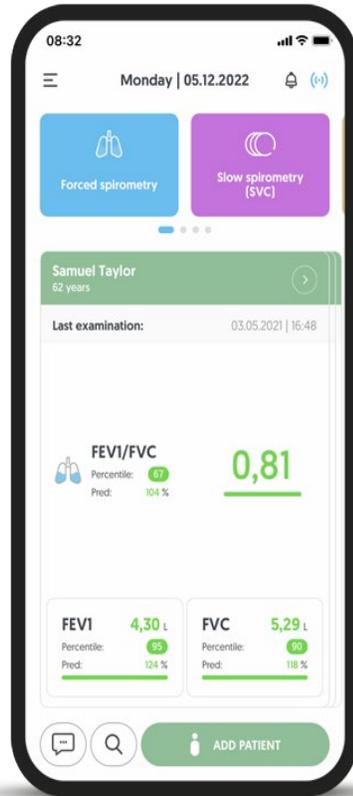

Creation of patients in the App for complete mobile testing

Pre and post bronchodilator workflow with prompt for administration of bronchodilator and selectable countdown timer for performing post test

Ability to perform SVC, FVC, PEF and SP02 measurements

Results instantly displayed along with F/V and V/T graphs and Z scores

Results can be securely exported in PDF format via email

Ability to add examination note by the healthcare professional overseeing the test



On the fly updates and new features - No Technician needed for any update.



Unlimited secure data storage and backup - Worry free and device independent operation.




Hassle free setup - Just runs in a browser and on a mobile platform.



Default secure and encrypted data transfer - always compliant.

# Software Panel

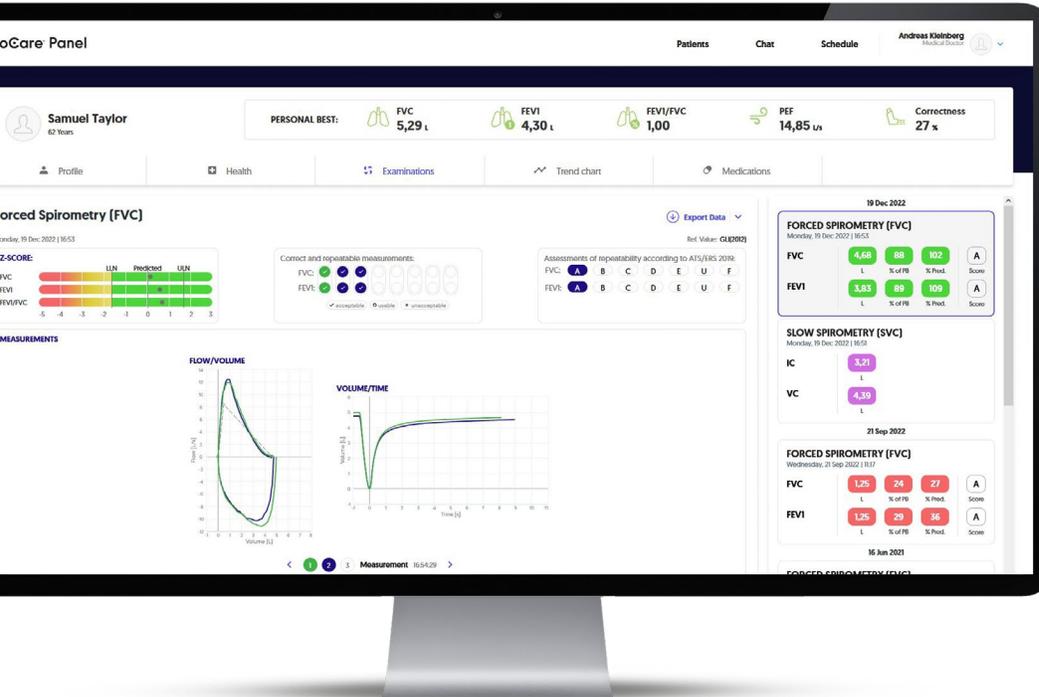


Instant access to all patient results on various platforms, irrespective of geographical location

Synchronization with App so results are instantly displayed

Ability to generate customised reports with Institute details on header. PDF report generation with reporting history and Export of results as XML, CSV or GDT

Clear concise quality grading for technically acceptable and reproducible maneuvers as per ATS/ERS 2019 Guidelines

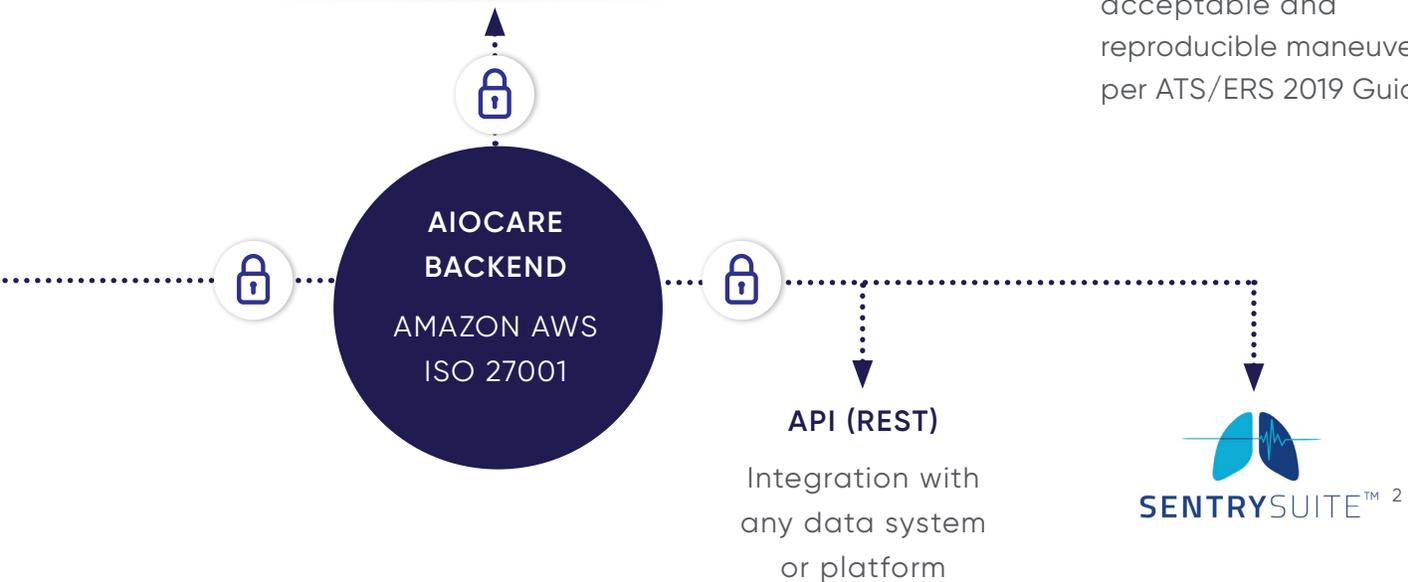

## REFERENCES

1. Based on the Bio Burden DIN EN ISO 11737-1: Report 18AA0088
2. Requires SentrySuite™ 3.30 or higher

The contents of this publication may differ from the current approval of the product or service in your country.

### GLOBAL HEADQUARTERS

Vyaire Medical, Inc.  
26125 N. Riverwoods Blvd.  
Mettawa, IL 60045  
USA

### AioCare

 HealthUp S. A.  
Twarda 18  
00-375 Warsaw:  
Poland

 2274

### MICROGARD™ II FILTER SENTRYSUITE™ SOFTWARE

 Vyaire Medical GmbH  
Leibnizstrasse 7  
97204 Hoechberg  
Germany

 0123

#### For international use only.

© 2023 Vyaire Medical, Inc. or one of its affiliates. All rights reserved. Vyaire and the Vyaire Medical logo are trademarks or registered trademarks of Vyaire Medical, Inc. or one of its affiliates. Please read the complete Instructions For Use that come with the devices or follow the instructions on the product labelling. Trademarks are the property of their respective owners VYR-INT-2000236 | 3.0

