



AioCare Doctor Better Diagnostics

Respiratory Disease Management



AioCare Modern Respiratory Disease Management.



Time efficient spirometry and diagnosis for experts and general practitioners alike.



Holistic patient treatment with smart home monitoring.



Wireless and truly mobile.



Decisions support algorithm for physicians and an AI-based quality assuring algorithm.



AioCare x MicroGard™ II Filter - Hygienic and secure.

"AioCare, an intelligent respiratory disease management to diagnose and monitor pulmonary diseases, even from the homes of my patients."





Hardware

Lung diagnostics & environment monitoring

AioCare

AioCare is a professional system to monitor and treat pulmonary diseases, consisting of a portable spirometer connected to a smartphone app and an online panel to access the results. The device enables a full spirometry test to be conducted in a doctor's office, including a bronchodilator test. Thanks to the applied algorithms, the results are instantly available on the phone and in a computer panel.


The smartphone app module supports the interpretation of the measured data. The system compares the results of subsequent spirometry tests before and after the application of the medicine, and automatically analyses the level and reversibility of obstruction diagnosis within 25 minutes.

AioCare can be used by patients at home to monitor their lungs' parameters, while doctors can track the results and consult if needed. Protecting their patients against the unexpected consequences of disease exacerbation.



Software

Mobile app & secure cloud-based data platform

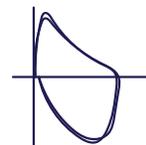

Data

Decisions support & AI-based quality assuring algorithms

Why AioCare?



Confirmed quality

A cutting edge, class IIa medical device (MDD), certified by the notified body TUV Nord. Validated by an observational study in a sample of 80 patients.



Accessibility

A small, portable device that can be taken anywhere. Safety tested by UL. Easy to use and inexpensive. Intuitive interface for creating patient records automatically saves test results, enables statistical analysis and management of the disease history.



4 easy steps to your first results


High accuracy of test results identical to those conducted with the medical equipment available at medical centres.



Warranty and support

2-year warranty and customer support line.

4 easy steps to your first results.

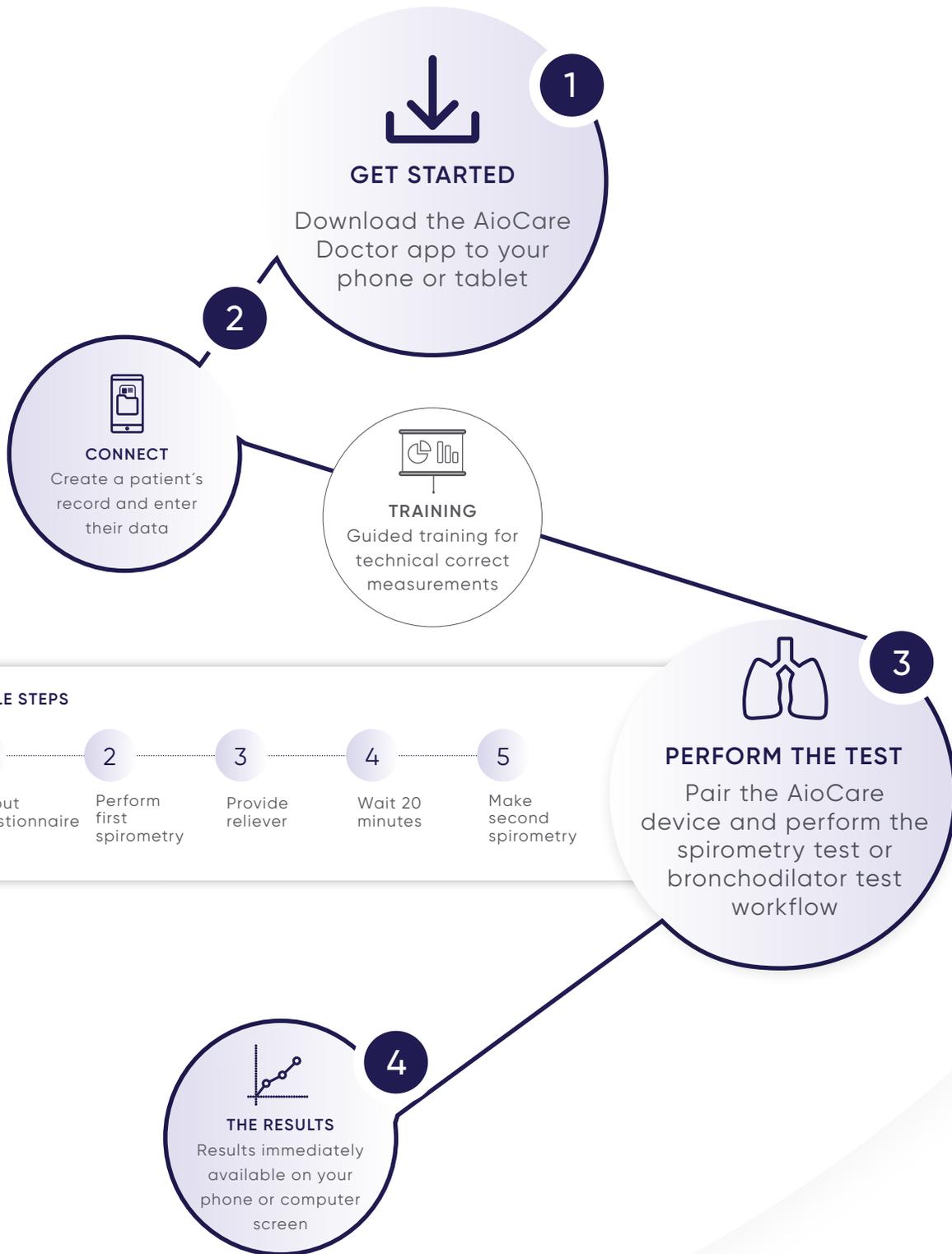

Full Spirometry

AioCare registers inspiratory & expiratory flow volume curves corrected by built-in real-time BTPS & displays key parameters, as well as SpO₂/HR.


Calibration free

Be ready for spirometry at any time.



CONNECTION BLUETOOTH 4.0


Connection secured by AES-128 encryption with Electronic CodeBook (ECB) and directly in the application by authorization of connection with the chosen AioCare spirometer.

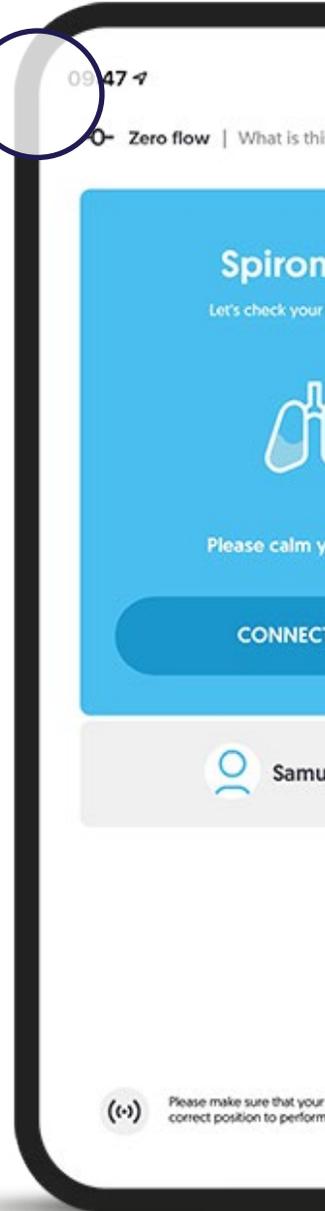

Simple Solution

Seamless IT integration.

Wireless connection & charging.

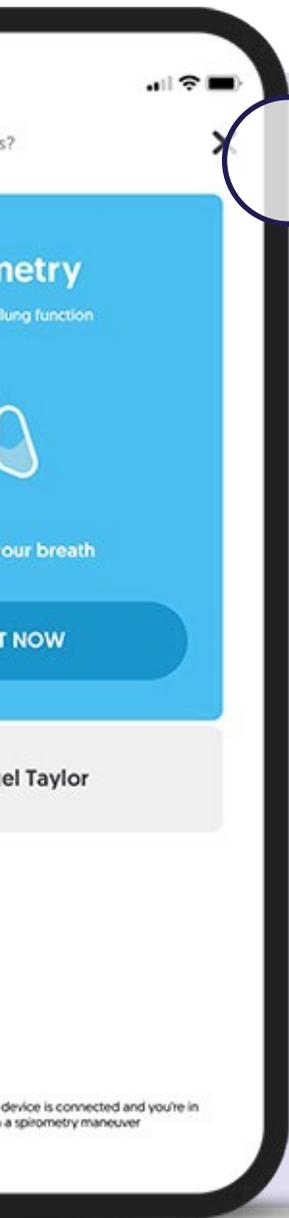

DATA TRANSFER
Transfer secured by SSL certificate (Secure Socket Layer)



API (REST)

Integration with any data system or platform



HIS (HL7)

Integration with Hospital Information System



SentrySuite® SOFTWARE

Profit from pulmonary expert reviews due to the peer review functionality in **SentrySuite**.
Enter full trend reports and combine with other pulmonary function test measures, as bodyplethysmography



AioCare on-line panel

Data aggregate in the desktop version for doctors



Better Diagnostics. Highest quality standards.



FEVI
Accuracy



FVC
Accuracy

ATS/ERS

Implemented 2019 guidelines



Patient Validation Tests



Multicenter Observation
Study



National obstruction
Diagnostic Program



Accuracy Validation
Report, May 2020

AioCare x MicroGard™ II Filter

High level validated cleaning methods

Protect patients, staff and device from cross-contamination by using the MicroGard II in-line filter.

MicroGard II filters provide 99.999% viral and bacterial efficiency against cross contamination (Nelson Test Report 1003754).



AioCare. Better lives.

AioCare means higher quality care for your patients. You give your patients a device to self-monitor and take care of their own health in a safe way, under your care.

Patients know their data is automatically shared with their doctor. It saves hours of telephone consultations, urgent visits and avoidable emergency interventions.

Stay connected with you patients, independent of any environmental influences.

An electronic journal enables the patient to avoid filling out tables by hand. The results of your AioCare patients are instantly visible in your panel. Additionally, you can easily access information about their medication and other symptoms. As your patients are no technicians, AioCare supports the user with additional sensors and a virtual guide to perform spirometry measurements in a technically correct way.




Disease control

individual digital asthma
self-control plan.




AioCare Respiratory Disease Management System

Doctor app
Transfer secured
by SSL certificate
(Secure Socket Layer)

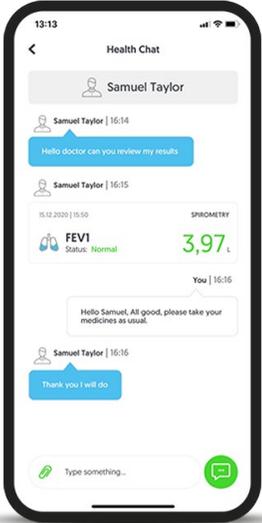

Patient app
Transfer secured
by SSL certificate
(Secure Socket Layer)

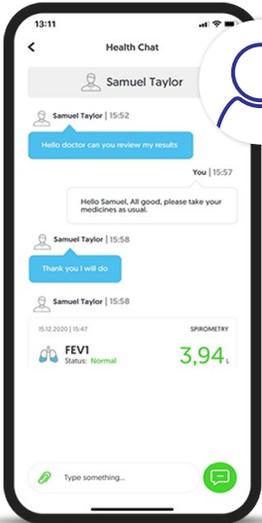


Contact initiated
by patient.
Measurement
results sent



Answer sent
by Doctor

Vyaire AioCare Services

Access to Professional Web Panel and Advanced App Features

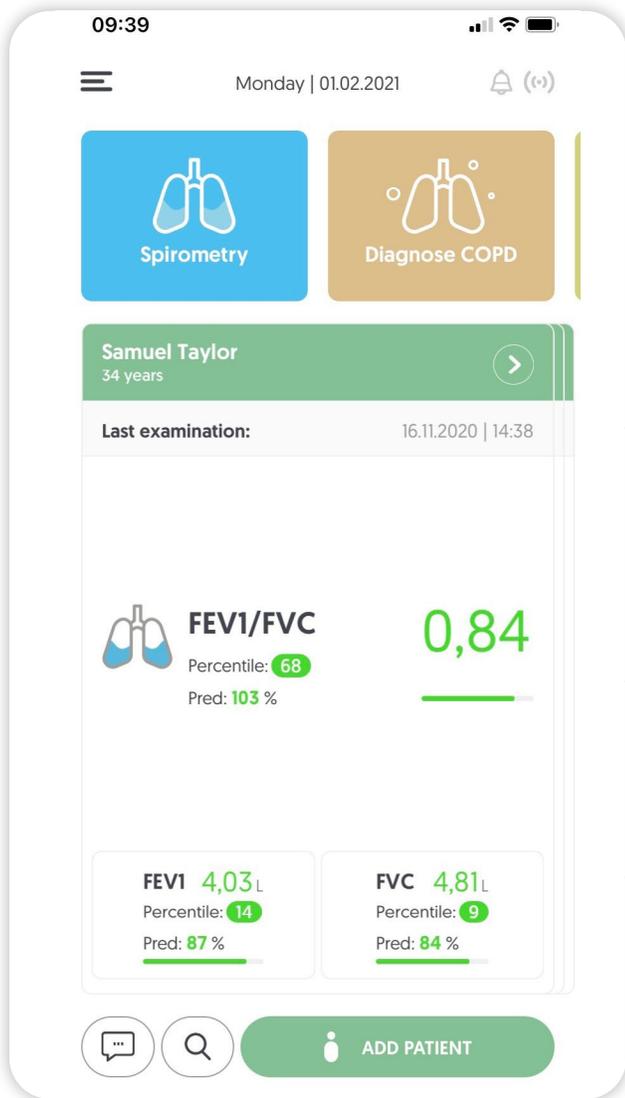


AioCare APP

Bronchodilator test with automatic data interpretation COPD diagnostic module

Advanced app feature: Access to patient home monitoring data

Export of results to PDF or Excel / share with email

Chat module for convenient 24/7 contact with your patient

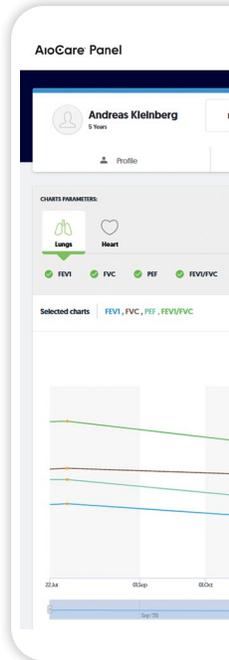

APP + PANEL:

Quick online and phone support of a dedicated technical specialist.



AioCare PANEL


Instant access to all results in an online web panel

Preview of data trends

Quick identification of exacerbations

Asthma Action Plan

Chat module for convenient 24/7 contact with your patient

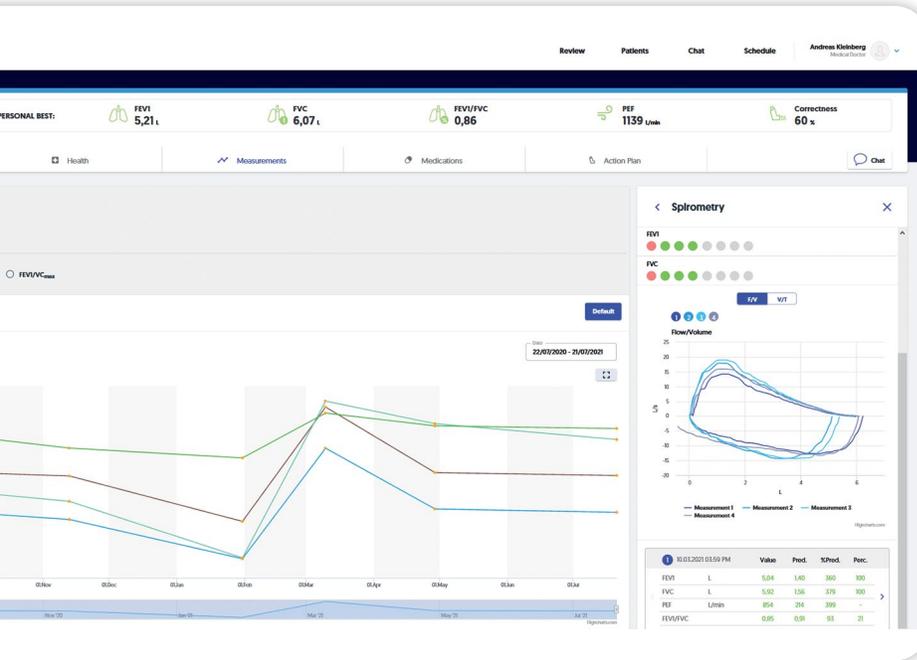

On the fly updates and new features - No Technician needed for any update.



Unlimited secure data storage and backup - Worry free and device independent operation.



Hassle free setup - Just runs in a browser and on a mobile platform.



Default secure and encrypted data transfer - always compliant.

REFERENCES

1. Based on the Bio Burden DIN EN ISO 11737-1: Report 18AA0088

GLOBAL HEADQUARTERS

Vyaire Medical, Inc.
26125 N. Riverwoods Blvd.
Mettawa, IL 60045
USA

AioCare

 HealthUp S. A.
Twarda 18
00-375 Warsaw:
Poland

 2274

MICROGARD™ II FILTER SENTRYSUITE® SOFTWARE

 Vyaire Medical GmbH
Leibnizstrasse 7
97204 Hoechberg
Germany

 0123

For Belgium, France, Germany, Italy, The Netherlands, Spain, Switzerland, United Kingdom, Spain, Saudi Arabia, Bahrain, UAE, Jordan, Kuwait, Oman distribution only.

© 2021 Vyaire Medical, Inc. or one of its affiliates. All rights reserved. Vyaire and the Vyaire Medical logo are trademarks or registered trademarks of Vyaire Medical, Inc. or one of its affiliates. Please read the complete Instructions For Use that come with the devices or follow the instructions on the product labelling. VYR-INT-2000236 v.2.0

© 2021 HealthUp SA. AioCare with logo is a registered trademark of HealthUp SA. For further details of HealthUp SA IP rights we refer to <https://aiocare.com/ip>.

